

A Cam-Med Company

### Brooks Stair Lift Proposal

arranty anlike manner and that 90 Days from the date atteries, any damage excludes any weather
anlike manner and that 90 Days from the date atteries, any damage excludes any weather
atteries, any damage excludes any weather
excludes any weather
\$465.00
\$450.00
\$420.00

# Make the most of your home with a Brooks Stairlift

Mr Frederick Brooks, a Daimler engineer, made his first stairlift in 1973. Almost 40 years later, by combining today's advanced technology with a tradition of reliability and comfort, we guarantee that the Brooks Lincoln Stairlift will help you make the most of your home.

Safety and comfort come first with the Brooks Stairlift. Self-charging DC power packs mean the lift is whisper quiet, incredibly smooth and will function even during a power cut. The slimline motor provides years of smooth, steady service, while its host of innovative features offer exceptional safety and comfort.

For over thirty years, we have provided the finest stairlifts to our customers. Ours is a heritage based on quality and service – a heritage you can trust. Join the thousands worldwide, who associate 'Brooks' with excellence.

We remain true to our guiding principles of quality and care. Here are just a few of the features that set us apart.

- Ultra quiet operation
- · Option to sit or stand
- Numerous Safety sensors
- Weatherproof outdoor stairlift available
- UL listed
- Register online for our comprehensive
   12-month warranty
- Order online 24/7
- Rapid response and delivery

Safety and comfort come first at Brooks. The stairlift uses DC power packs that are automatically recharged when the stairlift is at rest. This means that the lift is ultra quiet, incredibly smooth and will function even during a power outage.

"Ordering and receiving our Brooks Stairlift was far easier and quicker than I expected. Installation was clean and trouble-free. It's all been so simple and it has really given us back our freedom."







## Anatomy of the Brooks Stairlift


What makes the Brooks Lincoln the leading stairlift in the world? It is a combination of thoughtful design, top class engineering, reliability and user friendly features.

Every part of the stairlift has been designed with you in mind. Each feature is the result of careful consideration in response to the requirements of the stairlift user. We listen, and then we produce the best solution.

Brooks Stairlifts have slimline motors that provide a steady start / stop action and many years of faithful service. When you purchase a Brooks Lincoln Stairlift, you are buying freedom –

not costly maintenance. We are only a phone call away and our experts can give you all the advice you need.


Every member of the Brooks team has been trained to the highest standard. From your Authorized Brooks Dealer to the technician who services your stairlift, our care and attention to detail is second to none.

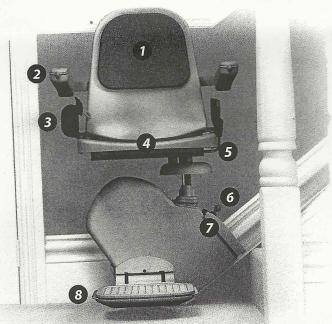

- 1 Padded seat and backrest.
- 2 Easy to use controls.
- 3 Secure safety belt.
- Folding arms, footrest and seat allow easy access to the stairway.
- Swivel seat no need for twisting the body to get on and off the lift.
- 6 Lockable on/off switch allows the user to prevent others from using the lift.
- Diagnostic digital display informs the user of the status of the lift.
- Five safety sensors stop the lift automatically if any obstruction is encountered.



Brooks, in keeping with its policy of continual development, reserves the right to change specifications without notice. All measurements are approximate. Images are shown for illustrative purposes only.

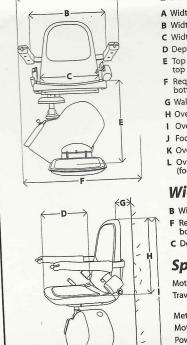

#### Dimensions

A Width outside armrest 22" 560mm
B Width inside armrest 17.5" 445mm
C Width of backrest 14.5" 370mm
D Depth of seat 16" 410mm

E Top of footrest to top of seat F Required space from bottom of stair

17.5" 445mm
26" 660mm

G Wall to seat width 1.5" 40mm

H Overall height of seat 19.5" 495mm

L Overall height 40" 1020mm

J Footrest surface to floor 4.5" 115mm
K Overall depth (folded) 13" 325mm

L Overall depth (footrest down)

#### **Wide Seat Option**

B Width inside armrest 20" 520mm
F Required space from bottom of stair
C Depth of seat 15" 380mm

#### Specifications

Motor power 0.5 hp
Travel speed No greater than 0.15 m/s

Method of drive Rack and Pinion Motor output torque 108 Nm

Power supply 24v DC (battery)
Maximum capacity 294 lbs (optional 350lbs)

Track Extruded Aluminium
Mains supply 110v, 1.5amp Charger

Dealer Stamp

Stairlift
gives the
user the
choice
whether to
sit or stand.

Sit/Stand



Perch Stairlift is ideal for customers who may have trouble sitting or bending at the knee.







